

Colony Counter

Two-Year Warranty

Features

Four digital displays

4 digit display offers convenience of showing the present and previous three count values.

Multi-counting average

Galaxy 330 can make average count from multiple counting.

Exclusive BACK button

The BACK button offers the convenience of deducting data when over counting.

Suitable for Petri-dishes in various diameters


Galaxy 330's adjustable dish holder can accept most Petri-dishes in a wide range of diameters from 90 to 150 mm or 60 mm.

Pressure sensor design

Pressure sensor design can allow user to operate with any marking pens. No worry about sample spoiled or contaminated.

Safety compliance

CE certification

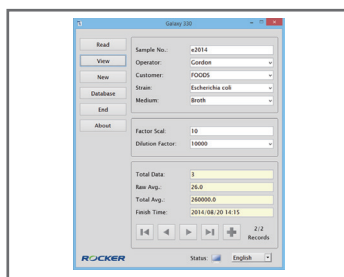

Galaxy 330, PC Software (optional)

In addition to recording and tracing data, the newest PC software can help user add, delete, search, print out and save data as Excel format (csv file).

Applications

- Total bacteria counting
- Total coliform counting

Ordering Information

175330-01

Galaxy 330, Colony Counter with U.S. plug (100-240V)

175330-02

Galaxy 330, Colony Counter with EU plug (100-240V)

175331-01

Galaxy 330, Colony Counter with software, U.S. plug (100-240V)

175331-02

Galaxy 330, Colony Counter with software, EU. plug (100-240V)

175330-12

Software Package for Galaxy330 (Include 8G USB flash and connecting cable)

Specification

Model	Galaxy 330
Display	3-digit (0-999) LED display, 4 panels
Magnification	2~3 X
Petri dish diameter	60 mm, 90-150 mm, adjustable
Memory	100
Average	Yes
Zero reset	Yes
Count back	Yes
Illumination	LED (adjustable)
Pressure sensitivity	Adjustable
Acoustical counting signal	Adjustable
Output	USB
Dimension (LxWxH)	35.1 x 30.2 x 37 cm
Net weight	5.6 kg